


## DUO-JET II

*Unlike traditional double lumen catheters, the Duo-Jet II utilizes a unique procedure first introduced by Canaud and later Tesio to allow exact placement of the cannulas in the vein before tunneling. This back-tunneling is facilitated by formable, steel trocars which are included in the sets.*

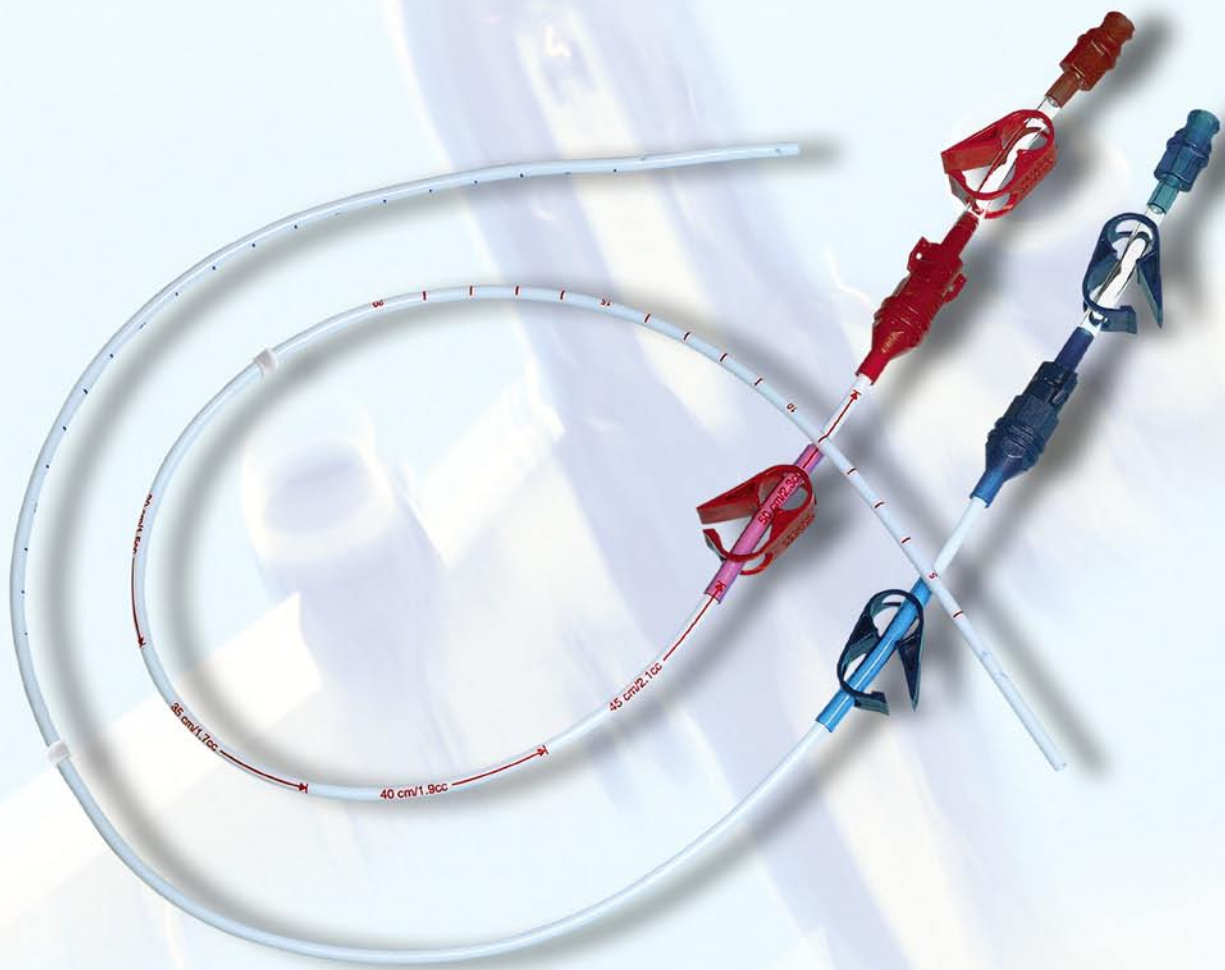

# DUO-JET II

## *features & benefits*

- Two-part adapters with extensions designate venous and arterial ports while offering ease of blood access.
- Two single lumen catheters allow blood flows of 400-450 ml/min.
- Two free-floating lumens independent of each other in the same vessel help minimize catheter occlusion by the vessel.
- Modified Seldinger approach provides easier and less traumatic access than surgical cut-down techniques.

SETS 5/BOX	DESCRIPTION
DJLT2000R	10FR x 52cm Duo-Jet II Catheter (right side approach) 22cm arterial, 25cm venous cuff placement
DJLT2000L	10FR x 52cm Duo-Jet II Catheter (left side approach) 27cm arterial, 30cm venous cuff placement

### SET CONTENTS:

- |                              |                        |
|------------------------------|------------------------|
| (2) Catheters                | (1) 6FR Sheath Dilator |
| (1) Arterial Extension (red) | (2) J/Flex Guidewires  |
| (1) Venous Extension (blue)  | (2) Injection Caps     |
| (2) Trocars                  | (1) StatLock®          |
| (2) Sheath Introducers       | (2) Clamps             |
| (2) Introducer Needles       |                        |

**Jet**  
Medical<sup>SA</sup>

Rue Fritz-Courvoisier 40  
2300 La Chaux-de-Fonds  
Switzerland  
Tel: +41-32-9679021  
Fax: +41-32-9679001  
www.jetmedical.net  
customerservice@jetmedical.net