

JETCATH[®]

- > HEMODIALYSIS
- > PLASMAPHERESIS
- > CONTINUOUS REPLACEMENT THERAPY



Jet
Medical[®] SA


JET CATH® APPLICATIONS

BRIDGE ACCESS - Awaiting creation/maturation of a permanent access with a graft.

TEMPORARY ACCESS - Waiting for a revision of a temporary malfunctioning fistula or graft.

ACUTE RENAL FAILURE - For renal replacement therapy in patients with multiple organ failure.

PLASMAPHERESIS - For immediate blood access.

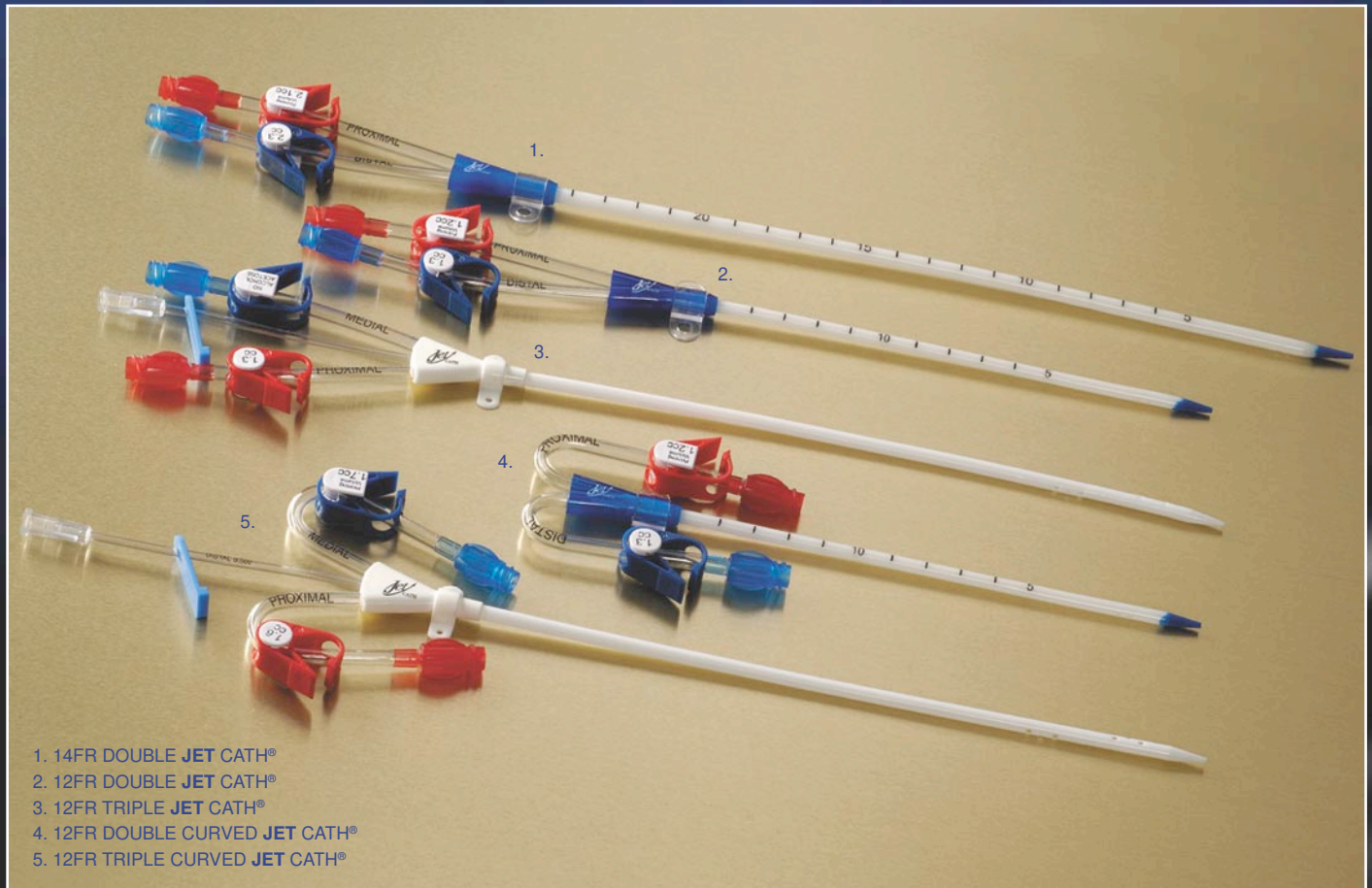

JET CATH® ZOOM IN > Take a closer look at the versatile Jet Cath® .



BIOCOMPATIBLE POLYURETHANE - Hemocompatible and resistant to breakdown caused by site care ointments and solutions. Tested for toxicity.

SOFT TIP - Reduces trauma to vessel wall upon introduction and in-situ. The radiopaque catheter and tip facilitates vessel location when visualized by X-ray.



TRIPLE LUMEN - Two dedicated arterial/venous access lumens and a separate central lumen for medications assures the integrity of the access lumens.

CURVED EXTENSIONS - The extensions are curved for jugular vein approaches to make the catheter more comfortable for the patient.



SOFT HUB - The soft hub allows for added patient comfort.

COLOR CODED LUERS - Venous and arterial luers color designed to resist cracking and breakdown over continued use.



FLOW RATES - The full line of "Jet Cath®" catheters provide superior flows with low pressures, maximizing treatment efficiency and effectiveness.

ORDERING INFORMATION

SETS 10/BOX	DESCRIPTION
DOUBLE STRAIGHT	
ARD1213S	13CM X 12FR Double Lumen Catheter
ARD1215S	15CM X 12FR Double Lumen Catheter
ARD1220S	20CM X 12FR Double Lumen Catheter
ARD1224S	24CM X 12FR Double Lumen Catheter
ARD1413S	13CM X 14FR Double Lumen Catheter
ARD1415S	15CM X 14FR Double Lumen Catheter
ARD1420S	20CM X 14FR Double Lumen Catheter
ARD1424S	24CM X 14FR Double Lumen Catheter
DOUBLE CURVED	
ARD1213C	13CM X 12FR Double Lumen Curved Catheter
ARD1215C	15CM X 12FR Double Lumen Curved Catheter
ARD1220C	20CM X 12FR Double Lumen Curved Catheter
ARD1224C	24CM X 12FR Double Lumen Curved Catheter
ARD1413C	13CM X 14FR Double Lumen Curved Catheter
ARD1415C	15CM X 14FR Double Lumen Curved Catheter
ARD1420C	20CM X 14FR Double Lumen Curved Catheter
ARD1424C	24CM X 14FR Double Lumen Curved Catheter
TRIPLE STRAIGHT	
ART1213S	13CM X 12FR Triple Lumen Catheter
ART1215S	15CM X 12FR Triple Lumen Catheter
ART1220S	20CM X 12FR Triple Lumen Catheter
ART1224S	24CM X 12FR Triple Lumen Catheter
TRIPLE CURVED	
ART1213C	13CM X 12FR Triple Lumen Curved Catheter
ART1215C	15CM X 12FR Triple Lumen Curved Catheter
ART1220C	20CM X 12FR Triple Lumen Curved Catheter
ART1224C	24CM X 12FR Triple Lumen Curved Catheter

DOUBLE LUMEN SETS CONTAIN:

- (1) CATHETER
- (1) INTRODUCER NEEDLE
- (1) VESSEL DILATOR
- (1) J-TIP GUIDEWIRE
- (2) INJECTION CAPS

TRIPLE LUMEN SETS CONTAIN:

- (1) CATHETER
- (1) INTRODUCER NEEDLE
- (1) VESSEL DILATOR
- (1) J-TIP GUIDEWIRE
- (3) INJECTION CAPS

Jet
Medical[®] SA

Rue Fritz-Courvoisier40
2300 La Chaux-de-Fonds
Switzerland
Tel: +41-32-9679021
Fax: +41-32-*9679001
www.jetmedical.net
customerservice@jetmedical.net