



AndraSnare® / LoopMaster Sochman Snare®

Gold plated atraumatic 3D nitinol snare /
Atraumatic 3D nitinol snare for foreign body retrieval

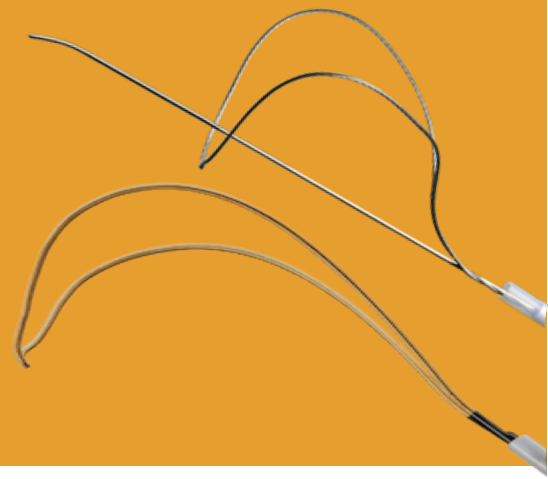

AndraSnare® & AndraSnare Micro®


1. Atraumatic nitinol loop with gold coating for better visualization in the X-rays
2. Pre-formed tip for better navigation
3. X-ray marker at the shaft output
4. Shaft visible in X-rays
5. Bend protection for shaft reinforcement at the Luer
6. Luer connection
7. Insertion part for simple insertion of the loop into the catheter
8. Large torquer for good operation
9. Continuous nitinol wire for 1:1 control of the loop

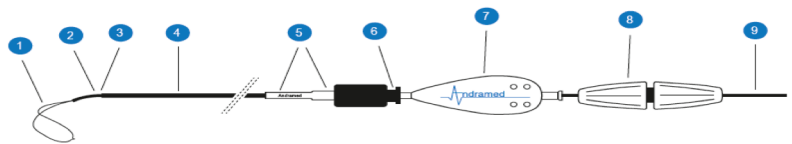

Order Codes AndraSnare® & AndraSnare Micro®

Order Codes	Loop Ø in mm	Vessel Ø in mm	Loop/Guidewire Length in cm	Introducer Sheath Length in cm	Introducer Sheath Ø in Fr
AS-5	5	3-5	125	110	4
AS-10	10	6-10	125	110	4
AS-15	15	11-15	125	110	4
AS-20	20	16-20	125	110	4
AS-25	25	21-25	125	110	5
AS-25/65	25	21-25	65	50	5
AS-25/150	25	21-25	150	135	5
AS-30	30	26-30	125	110	5
AS-35	35	30-35	125	110	5
ASM-2	2	1-2	175	150	2,3-3
ASM-4	4	3-4	175	150	2,3-3
ASM-7	7	5-7	175	150	2,3-3

LoopMaster® Sochman Snare

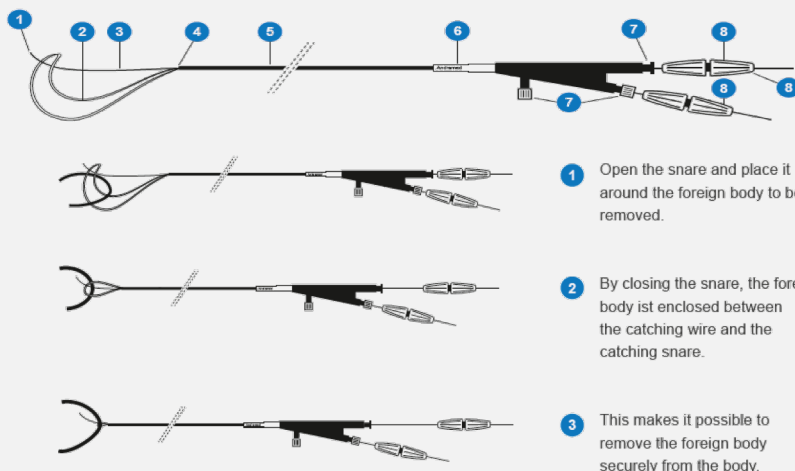

1. Atraumatic distal end at the catching wire
2. Nitinol loop
3. Catching wire
4. X-ray marker
5. Flexible shaft
6. Bend protection at the Luer
7. Luer connector
8. Large torquers for simple control of the catching snare and the catching wire

Order Code LoopMaster® Sochman Snare:
ASL-25